

Customer and Technical Service: 800-822-2947

July 2006

PN: 400-7137 Rev.: D

© 2003, Abaxis, Inc., Union City, CA 94587

## 1. Intended Use

The Piccolo<sup>®</sup> Hepatic Function Panel reagent disc, used with the Piccolo Blood Chemistry Analyzer, is intended to be used for the *in vitro* quantitative determinations of alanine aminotransferase, albumin, alkaline phosphatase, aspartate aminotransferase, direct bilirubin, total bilirubin, and total protein in heparinized whole blood, heparinized plasma, or serum in a clinical laboratory setting.

## 2. Summary and Explanation of Tests

The Piccolo Hepatic Function Panel reagent disc and the Piccolo Blood Chemistry Analyzer comprise an *in vitro* diagnostic system that aids the physician in diagnosing the following disorders.

Alanine aminotransferase:	Liver diseases, including viral hepatitis and cirrhosis; heart diseases.
Albumin:	Liver and kidney diseases.
Alkaline phosphatase:	Liver, bone, parathyroid, and intestinal diseases.
Aspartate aminotransferase:	Liver disease including hepatitis and viral jaundice, shock.
Direct bilirubin:	Liver disorders, hemolytic hematological, and metabolic disorder, including hepatitis and gall bladder obstruction.
Total bilirubin:	Liver disorders, including hepatitis and gall bladder obstruction; jaundice.
Total protein:	Liver, kidney, bone marrow diseases; metabolic and nutritional disorders.


**As with any diagnostic test procedure, all other test procedures including the clinical status of the patient should be considered prior to final diagnosis.**

## 3. Principle of Procedure

### Alanine Aminotransferase (ALT)

Alanine aminotransferase (ALT) has been measured by three methods. Two of these methods—the colorimetric dinitrophenylhydrazine coupling technique<sup>1,2</sup> and the fluorescent enzymatic assay—are rarely used.<sup>3</sup> An enzymatic method based on the work of Wróblewski and LaDue<sup>4</sup> is the most common technique for determining ALT concentrations in serum. A modified Wróblewski and LaDue procedure has been proposed as the recommended procedure of the International Federation of Clinical Chemistry (IFCC).<sup>5</sup>

The method developed for use on the Piccolo Analyzer is a modification of the IFCC-recommended procedure. In this reaction, ALT catalyzes the transfer of an amino group from L-alanine to  $\alpha$ -ketoglutarate to form L-glutamate and pyruvate. Lactate dehydrogenase catalyzes the conversion of pyruvate to lactate. Concomitantly, NADH is oxidized to NAD<sup>+</sup>, as illustrated in the following reaction scheme.

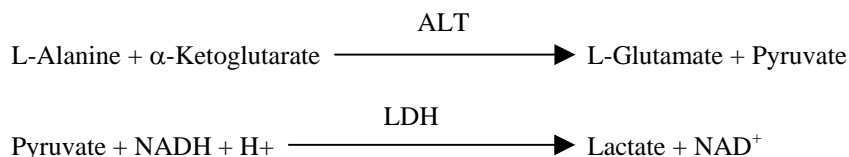


The rate of change of the absorbance difference between 340 nm and 405 nm is due to the conversion of NADH to NAD<sup>+</sup> and is directly proportional to the amount of ALT present in the sample.

### Albumin (ALB)

Early methods used to measure albumin include fractionation techniques<sup>6-8</sup> and tryptophan content of globulins.<sup>9,10</sup> These methods are unwieldy to perform and do not have a high specificity. Two immunochemical techniques are considered as reference methods, but are expensive and time consuming.<sup>11</sup> Dye binding techniques are the most frequently used methods for measuring albumin. Bromocresol green (BCG) is the most commonly used of the dye binding methods but may over-estimate albumin concentration, especially at the low end of the normal range.<sup>12</sup> Bromocresol purple (BCP) is the most specific of the dyes in use.<sup>13,14</sup>

Bromocresol purple (BCP), when bound with albumin, changes color from a yellow to blue color. The absorbance maximum changes with the color shift.

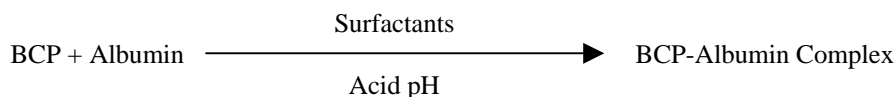


Bound albumin is proportional to the concentration of albumin in the sample. This is an endpoint reaction that is measured as the difference in absorbance between 600 nm and 550 nm.

### Alkaline Phosphatase (ALP)

Techniques to measure alkaline phosphatase were first developed over 60 years ago. Several of these endpoint or two-point spectrophotometric methods<sup>15,16</sup> are now considered obsolete or too cumbersome. The use of *p*-nitrophenyl phosphate (*p*-NPP) increased the speed of the reaction.<sup>17,18</sup> The reliability of this technique was greatly increased by the use of a metal-ion buffer to maintain the concentration of magnesium and zinc ions in the reaction.<sup>19</sup> The American Association for Clinical Chemistry (AACC) reference method<sup>20</sup> uses *p*-NPP as a substrate and a metal-ion buffer.

The Piccolo procedure is modified from the AACC and IFCC<sup>21</sup> methods. Alkaline phosphatase hydrolyzes *p*-NPP in a metal-ion buffer and forms *p*-nitrophenol and phosphate.

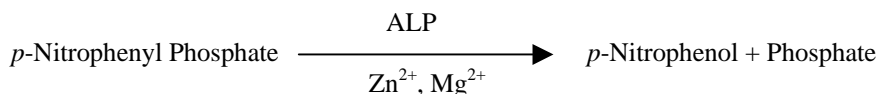


The amount of ALP in the sample is proportional to the rate of increase in absorbance difference between 405 nm and 500 nm.

### Aspartate Aminotransferase (AST)

The aspartate aminotransferase (AST) test is based on the Karmen rate method<sup>22</sup> as modified by Bergmeyer.<sup>23</sup> The current International Federation of Clinical Chemistry (IFCC) reference method utilizes the Karmen / Bergmeyer technique of coupling malate dehydrogenase (MDH) and reduced nicotinamide dinucleotide (NADH) in the detection of AST in serum.<sup>23,24</sup> Lactate dehydrogenase (LDH) is added to the reaction to decrease interference caused by endogenous pyruvate.

AST catalyzes the reaction of L-aspartate and  $\alpha$ -ketoglutarate into oxaloacetate and L-glutamate. Oxaloacetate is converted to malate and NADH is oxidized to NAD<sup>+</sup> by the catalyst MDH.

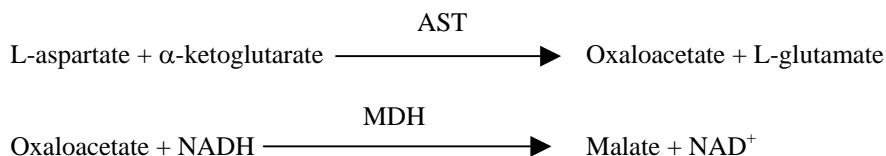


The rate of absorbance change at 340 nm / 405 nm caused by the conversion of NADH to NAD<sup>+</sup> is directly proportional to the amount of AST present in the sample.

### Direct Bilirubin (DBIL)

Direct Bilirubin was first detected in human serum by Van Den Bergh and Müller when they observed that the pigment in human bile reacted with diazotized sulfanilic acid in the absence of alcohol.<sup>25</sup> Today, commonly used techniques to measure direct bilirubin are modifications of a method developed by Malloy and Evelyn<sup>26</sup> which relies on diazotized sulfanilic acid combining

with bilirubin to form the chromophore azobilirubin. Some diazo methods produce unreliable results because unconjugated may be quantitated as direct bilirubin.<sup>27</sup> A more specific assay for total bilirubin was developed after the enzyme bilirubin oxidase was isolated from *Myrothecium verrucaria* MT-1.<sup>28-30</sup> This method is also specific for direct bilirubin when the reaction is run at a lower pH.<sup>31,32</sup>

In the enzymatic procedure, the soluble bilirubin complex (direct bilirubin) is oxidized by bilirubin oxidase into biliverdin.

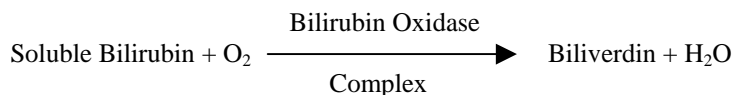


Direct Bilirubin is quantitated as the difference in absorbance between 467nm and 550nm. The initial absorbance of this end point reaction is determined from the direct bilirubin blank cuvette and the final absorbance is obtained from the direct bilirubin test cuvette. The amount of direct bilirubin in the sample is proportional to the difference between the initial and final absorbance measurements.

#### Total Bilirubin (TBIL)

Total bilirubin levels have been typically measured by tests that employ diazotized sulfanilic acid.<sup>26,33</sup> A newer, more specific method has been developed using the enzyme bilirubin oxidase.<sup>28-30</sup> In addition to using the more specific total bilirubin test method, photodegradation of the analyte is minimized on the Piccolo Analyzer because the sample can be tested immediately after collection.

In the enzyme procedure, bilirubin is oxidized by bilirubin oxidase into biliverdin.

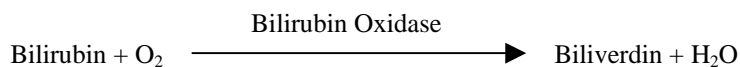


Bilirubin is quantitated as the difference in absorbance between 467 nm and 550 nm. The initial absorbance of this endpoint reaction is determined from the bilirubin blank cuvette and the final absorbance is obtained from the bilirubin test cuvette. The amount of bilirubin in the sample is proportional to the difference between the initial and final absorbance measurements.

#### Total Protein (TP)

The total protein method is a modification of the biuret reaction, noted for its precision, accuracy, and specificity.<sup>34</sup> Originally developed by Riegler<sup>35</sup> and modified by Weichselbaum<sup>36</sup>, Dumas, et al.<sup>37</sup> proposed a biuret reaction as a candidate total protein reference method.

In the biuret reaction, the protein solution is treated with cupric [Cu(II)] ions in a strong alkaline medium. Sodium potassium tartrate and potassium iodide are added to prevent the precipitation of copper hydroxide and the auto-reduction of copper, respectively.<sup>36</sup> The Cu(II) ions react with peptide bonds between the carbonyl oxygen and amide nitrogen atoms to form a colored Cu-protein complex.

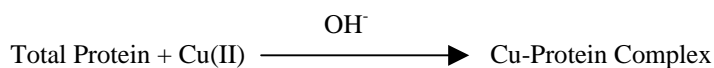

The amount of total protein present in the sample is directly proportional to the absorbance of the Cu-protein complex. The total protein test is an endpoint reaction and the absorbance is measured as the difference in absorbance between 550 nm and 850 nm.

### 4. Principle of Operation

See the Piccolo Blood Chemistry Analyzer Operator's Manual, for the Principles and Limitations of the Procedure.

### 5. Description of Reagents

#### Reagents

Each Piccolo Hepatic Function Panel reagent disc contains dry test-specific reagent beads (described below). A dry sample blank reagent (comprised of buffer, surfactants, excipients, and preservatives) is included in each disc for use in calculating concentrations of alanine aminotransferase (ALT), albumin (ALB), alkaline phosphatase (ALP), aspartate aminotransferase (AST). Dedicated sample blanks are included in the disc for total bilirubin (TBIL), direct bilirubin (DBIL), and total protein (TP). Each reagent disc also contains a diluent consisting of surfactants, excipients, and preservatives.

**Table 1: Reagents**

---

Component	Quantity/Disc
L-Alanine	874 µg
L-Aspartic acid	426 µg
Bilirubin oxidase	0.2 U
Bromcresol purple	2 µg
Cupric sulfate	134 µg
α-Ketoglutaric acid	82 µg
Lactate dehydrogenase	0.13 U
Magnesium chloride	3 µg
Malate dehydrogenase (MDH) (porcine heart)	0.01 U
β-Nicotinamide adenine dinucleotide, reduced (NADH)	12 µg
p-NPP	56 µg
Potassium iodide	28 µg
Sodium potassium tartrate	343 µg
Zinc sulfate	3 µg
Buffers, Surfactants, excipients, and preservatives	

**Warnings & Precautions**

- For *in vitro* Diagnostic Use
- The diluent container in the reagent disc is automatically opened when the analyzer drawer closes. A disc with an opened diluent container cannot be re-used. Ensure that the sample or control has been placed into the disc before closing the drawer
- Used reagent discs contain human body fluids. Follow good laboratory safety practices when handling and disposing used discs.<sup>38</sup> See the Piccolo Blood Chemistry Analyzer Operator's Manual for instructions on cleaning biohazardous spills.
- The reagent discs are plastic and may crack or chip if dropped. **Never** use a dropped disc as it may spray biohazardous material throughout the interior of the analyzer.
- Reagent beads may contain acids or caustic substances. The operator does not come into contact with the reagent beads when following the recommended procedures. In the event that the beads are handled (e.g., cleaning up after dropping and cracking a reagent disc), avoid ingestion, skin contact, or inhalation of the reagent beads.

**Instructions for Reagent Handling**

Reagent discs may be used directly from the refrigerator without warming. Do not allow discs sealed in foil pouches to remain at room temperature longer than 48 hours prior to use. Open the sealed foil pouch, remove the disc, and use according to the instructions provided in the Piccolo Blood Chemistry Analyzer Operator's Manual. A disc not used within 20 minutes of opening the pouch should be discarded.

**Storage**

Store reagent discs in their sealed pouches at 2-8°C (36-46°F). Do not expose opened or unopened discs to direct sunlight or temperatures above 32°C (90°F). Reagent discs may be used until the expiration date included on the package. The expiration date is also encoded in the bar code printed on the bar code ring. An error message will appear on the Piccolo Blood Chemistry Analyzer display if the reagents have expired.

**Indications of Reagent Disc Instability/Deterioration**

A torn or otherwise damaged pouch may allow moisture to reach the unused disc and adversely effect reagent performance. Do not use a rotor from a damage pouch.

**6. Instrument**

See the Piccolo Blood Chemistry Analyzer Operator's Manual for complete information on using the analyzer.

## 7. Sample Collection and Preparation

Sample collection techniques are described in the Piccolo Blood Chemistry Analyzer Operator's Manual.

- The minimum required sample size is ~100 µL of heparinized whole blood, heparinized plasma, serum, or serum control. The reagent disc sample chamber can contain up to 120 µL of sample.
- Whole blood samples obtained by venipuncture must be homogeneous before transferring a sample to the reagent disc. Gently invert the collection tube several times just prior to sample transfer. Do **not** shake the collection tube; shaking can cause hemolysis.
- Whole blood venipuncture samples should be run within 60 minutes of collection.<sup>39</sup> Refrigerating whole blood samples can cause significant changes in concentrations of **aspartate aminotransferase**.<sup>40</sup> The sample may be separated into plasma or serum and stored in capped sample tubes at 2-8°C (36-46°F) if the sample can not be run within 60 minutes.
- **Total and direct bilirubin** results may be adversely affected by photodegradation.<sup>41</sup> Whole blood samples not run immediately should be stored in the dark for no longer than 60 minutes. If the sample can not be analyzed within that period, it should be separated into plasma or serum and stored in a capped sample tube in the dark at low temperatures.<sup>42</sup>
- Use only lithium heparin evacuated specimen collection tubes for whole blood or plasma samples. Use no-additive evacuated specimen collection tubes or serum separator tubes for serum samples.

## 8. Procedure

### Materials Required

- One Piccolo Hepatic Function Panel Reagent Disc

### Materials Required but Not Provided

- Piccolo Blood Chemistry Analyzer
- Commercially available control reagents recommended by Abaxis (refer to the Piccolo Blood Chemistry Analyzer Operator's Manual).

### Test Parameters

- The Piccolo Blood Chemistry Analyzer operates at ambient temperatures between 15°C and 32°C (59-90°F). The analysis time for each Piccolo Hepatic Function Panel reagent disc is less than 14 minutes. The analyzer maintains the reagent disc at a temperature of 37°C (98.6°F) over the measurement interval.

### Test Procedure

The complete sample collection and step-by-step operating procedures are detailed in the Piccolo Blood Chemistry Analyzer Operator's Manual.

### Calibration

The Piccolo Blood Chemistry Analyzer is calibrated by the manufacturer before shipment. The bar code printed on the reagent disc bar code ring provides the analyzer with disc-specific calibration data. See the Piccolo Analyzer Operator's Manual.

### Quality Control

Performance of the Piccolo Blood Chemistry Analyzer can be verified by running controls. Controls recommended by Abaxis are listed in the Piccolo Blood Chemistry Analyzer Operator's Manual. Other human serum or plasma-based controls may not be compatible.

See the Piccolo Analyzer Operator's Manual, for a detailed discussion on running, recording, interpreting, and plotting control results.

## 9. Results

The Piccolo Blood Chemistry Analyzer automatically calculates and prints the analyte concentrations in the sample. Details of the endpoint and rate reaction calculations are found in the Piccolo Blood Chemistry Analyzer Operator's Manual.

Interpretation of results is detailed in the Piccolo Blood Chemistry Analyzer Operator's Manual. Results are printed onto results cards supplied by Abaxis. The result cards have an adhesive backing for easy placement in the patient's files.

## 10. Limitations of Procedure

General procedural limitations are discussed in the Piccolo Blood Chemistry Analyzer Operator's Manual.

- The only anticoagulant **recommended for use** with the Piccolo Blood Chemistry System is **lithium heparin**. Abaxis has performed studies demonstrating that EDTA, fluoride, oxalate, and any anticoagulant containing ammonium ions will interfere with at least one chemistry contained in the Piccolo Hepatic Function Panel Reagent Disc.
- Samples with hematocrits in excess of 62% packed red cell volume may give inaccurate results. Samples with high hematocrits may be reported as hemolyzed. These samples may be spun down to get plasma then re-run in a new reagent disc.
- **Any result for a particular test that exceeds the assay range should be analyzed by another approved test method or sent to a referral laboratory. Do not dilute the sample and run it again on the Piccolo Blood Chemistry Analyzer.**

**Warning:** Extensive testing of the Piccolo Blood Chemistry System has shown that in very rare instances, sample dispensed into the reagent disc may not flow smoothly into the sample chamber. Due to uneven flow, an inadequate quantity of sample may be analyzed and several results may fall outside the reference ranges. The sample may be re-run using a new reagent disc.

### Interference

Substances were tested as interferents with the analytes. Human serum pools were prepared. The concentration at which each potential interferent was tested was based on testing levels in NCCLS EP7-A.<sup>45</sup>

### Effects of Endogenous Substances

- Physiological interferents (hemolysis, lipemia, and icterus) cause changes in the reported concentrations of some analytes. The sample indices are printed on the bottom of each test card to inform the operator about the levels of interferents present in each sample. The Piccolo Blood Chemistry System suppresses any results that are affected by > 10 % interference from hemolysis, icterus, and lipemia. "HEM", "LIP", "ICT" respectively, is printed on the results card in place of the result.

### Effects of Therapeutic Substances

- Significant interference is defined as a >10% shift in the results for a normal range specimen. Human pools were supplemented with a known concentration of the drugs or chemicals and then analyzed.

**Table 2: Therapeutic Substances Evaluated**

	Physiologic or Therapeutic Range <sup>43-48</sup> (mg/dL)	Highest Concentration Tested (mg/dL)
<b>Endogenous Substances</b>		
Acetaminophen	1-2	100
Acetylsalicylic acid	2-10	50
Chloramphenicol	1-2.5	100
Cimetidine	0.1-1	16
Erythromycin	0.2-2	10
Isoniazide	0.1-0.7	4
Ketoprofen	—	50
Methicillin	—	100
Methotrexate	0.1	0.5
Metronidazole	0.1	5
Nafcillin	—	1
Oxacillin	—	1
Phenytoin	1-2	3

**Table 3: Substances With Significant Interference > 10%**

	<b>Physiologic or Therapeutic Range<sup>43-48</sup> (mg/dL)</b>	<b>Concentration with No Significant Interference (mg/dL)</b>	<b>Interference<sup>A</sup></b>
<b>Alanine Aminotransferase (ALT)</b>			
Ascorbic acid	0.8-1.2	20	11% inc
Oxaloacetate	—	132	843% inc
<b>Albumin (ALB)</b>			
Acetoacetate	0.05-3.60	102	18% dec
Ampicillin	0.5	30	12% dec
Caffeine	0.3-1.5	10	14% dec
Calcium chloride	—	20	17% dec
Cephalothin (Keflin)	10	400	13% inc
Ibuprofen	0.5-4.2	50	28% inc
α-Ketoglutarate	—	5	11% dec
Nitrofurantoin	0.2	20	13% dec
Proline	—	4	12% inc
Sulfalazine	2-4	10	14% dec
Sulfanilamide	10-15	50	12% dec
Theophylline	1-2	20	11% dec
<b>Alkaline Phosphatase (ALP)</b>			
Theophylline	1-2	20	42% dec
<b>Aspartate Aminotransferase (AST)</b>			
	None	None	None
<b>Direct Total Bilirubin (DBIL)</b>			
Ascorbic acid	0.8-1.2	2.5	30% dec
Dopamine	0.3-1.5	15	50% dec
<b>Total Bilirubin (TBIL)</b>			
Dopamine	—	19	55% dec
L-dopa	—	5	17% dec
<b>Total Protein (TP)</b>			
	None	None	None

<sup>A</sup> Dec.= decreased concentration of the specified analyte; Inc. = increased concentration of the specified analyte

For additional information on potential chemical interferents, see the Bibliography.

## 11. Expected Values

Samples from a total of 125 adult males and females analyzed on the Piccolo Blood Chemistry Analyzer, were used to determine the reference intervals. These ranges are provided as a guideline only. ALP levels in growing children are highly variable.<sup>47</sup> It is recommended that your office or institution establish normal ranges for your particular patient population.

**Table 4: Piccolo References Intervals**

Analyte	Reference Interval	
	Common Units	SI Units
Alanine Amino-transferase (ALT)	10-47 U/L	10-47 U/L
Albumin (ALB)	3.3-5.5 g/dL	33-55 g/L
Alkaline Phosphatase (ALP), Male	53-128 U/L	53-128 U/L
Alkaline Phosphatase (ALP), Female	42-141 U/L	42-141 U/L
Aspartate Amino-transferase (AST)	11-38 U/L	11-38 U/L
Direct Bilirubin (DBIL)	0-0.3 mg/dL	0-5.1 µmol/L
Total Bilirubin (TBIL)	0.2-1.6 mg/dL	3.4-27.4 µmol/L
Total Protein (TP)	6.4-8.1 g/dL	64-81 g/L

**12. Performance Characteristics****Linearity**

The chemistry for each analyte is linear over the dynamic range listed below when the Piccolo Blood Chemistry Analyzer is operated according to the recommended procedure (refer to the Piccolo Blood Chemistry Analyzer Manual).

**Table 5: Piccolo Dynamic Ranges**

Analyte	Dynamic Range	
	Common Units	SI Units
Alanine Aminotransferase (ALT)	5-2000 U/L	5-2000 U/L
Albumin (ALB)	1-6.5 g/dL	10-65 g/L
Alkaline Phosphatase (ALP)	5-2400 U/L	5-2400 U/L
Aspartate Amino-transferase (AST)	5-2000 U/L	5-2000 U/L
Direct Bilirubin (DBIL)	0.1-15 mg/dL	1.7-257 µmol/L
Total Bilirubin (TBIL)	0.1-30 mg/dL	1.7-513 µmol/L
Total Protein (TP)	2-14 g/dL	20-140 g/L

**Specificity**

The lower limit of detection for each analyte is: alanine aminotransferase 5 U/L; albumin 1 g/dL (10 g/L); alkaline phosphatase 5 U/L; aspartate aminotransferase 5 U/L; direct bilirubin 0.1mg/dL (1.7 µmol/L); total bilirubin 0.1 mg/dL (1.7 µmol/L); and total protein 2 g/dL (20 g/L).

**Precision**

For all assays precision studies were conducted using NCCLS EP5-A guidelines.<sup>49</sup> Results for within-run and total precision were determined by testing two levels of control material. Controls were run in duplicate twice each day for 20 days over a four-week period. Results of the precision studies are shown in Table 6.

**Table 6: Precision (N=80)**

Analyte	Within-Run	Total
<b>Alanine Aminotransferase (U/L)</b>		
<u>Control Level 1</u>	21	21
Mean	2.76	2.79
SD	13.4	13.5
%CV		
<u>Control Level 2</u>		
Mean	52	52
SD	2.70	3.25
%CV	5.2	6.2
<b>Albumin (g/dL)</b>		
<u>Control Level 1</u>		
Mean	5.6	5.6

**Table 6: Precision (N=80) (continued)**

Analyte	Within-Run	Total
SD	0.09	0.11
%CV	1.7	2.1
<u>Control Level 2</u>		
Mean	3.7	3.7
SD	0.07	0.11
%CV	2.0	2.9
<b>Alkaline Phosphatase (U/L)</b>		
<u>Control Level 1</u>		
Mean	39	39
SD	1.81	2.29
%CV	4.6	5.8
<u>Control Level 2</u>		
Mean	281	281
SD	4.08	8.75
%CV	1.5	3.1
<b>Aspartate Aminotransferase (U/L)</b>		
<u>Control Level 1</u>		
Mean	47	49
SD	0.98	0.92
%CV	2.1	1.9
<u>Control Level 2</u>		
Mean	145	147
SD	1.83	1.70
%CV	1.3	1.2
<b>Direct Bilirubin (mg/dL)</b>		
<u>Control Level 1</u>		
Mean	0.4	0.4
SD	0.03	0.03
%CV	6.5	6.6
<u>Control Level 2</u>		
Mean	2.2	2.2
SD	0.10	0.12
%CV	4.8	5.6
<b>Total Bilirubin (mg/dL)</b>		
<u>Control Level 1</u>		
Mean	0.8	0.8
SD	0.06	0.07
%CV	8.0	9.3
<u>Control Level 2</u>		
Mean	5.2	5.2
SD	0.09	0.15
%CV	1.7	2.8
<b>Total Protein (g/dL)</b>		
<u>Control Level 1</u>		
Mean	6.8	6.8
SD	0.05	0.08
%CV	0.8	1.2
<u>Control Level 2</u>		
Mean	4.7	4.7
SD	0.09	0.09
%CV	2.0	2.0

**Correlation**

Heparinized whole blood and serum samples were collected from patients at two sites. The whole blood samples were analyzed by the Piccolo Blood Chemistry Analyzer at the field sites and the serum samples were analyzed by the Piccolo Analyzer and by comparative methods. In some cases, high and low supplemented samples were used to cover the dynamic range. All samples were run in singlicate on the same day. Representative correlation statistics are shown in Table 7.

**Table 7: Correlation of Piccolo Blood Chemistry Analyzer with Comparative Method**

Analyte	Correlation coefficient	Slope	Intercept	SEE	N	Sample Range	Comparative Method
<b>Alanine Aminotransferase</b> (U/L)	0.981	0.905	1.3	3.21	86	10-174	Paramax®
	0.985	0.946	-2.5	2.84	67	10-174	Technicon
<b>Albumin</b> (g/dL)	0.854	1.001	-0.3	0.22	261	1.1-5.3	Paramax®
	0.896	0.877	-0.1	0.21	100	1.5-5.0	Beckman
<b>Alkaline Phosphatase</b> (U/L)	0.988	0.970	-5.9	3.97	99	27-368	Paramax®
	0.929	1.136	-17.6	4.79	80	26-150	Technicon
<b>AspartateAminotransferase</b> (U/L)	0.93	0.87	5.3	2.76	159	13-111	Paramax®
	1.0	0.97	3.0	1.90	46	13-252	DAX™
<b>Direct Bilirubin</b> (mg/dL)	0.990	0.88	-0.1	0.08	263	0-12.8	Paramax®
<b>Total Bilirubin</b> (mg/dL)	0.974	0.901	0.0	0.07	250	0.2-3.7	Paramax®
	0.980	1.113	-0.4	0.09	91	0.1-6.4	Beckman
<b>Total Protein</b> (g/dL)	0.849	0.932	0.6	0.19	251	5.7-9.2	Paramax®
	0.873	0.935	0.3	0.16	92	6.5-9.2	Beckman

\*Serum samples from hospitalized patients provided a broader, and possibly more useful, sample range than did venous whole blood samples from outpatients.

### 13. Bibliography

1. Tonhazy NE, White NG, Umbreit WW. A rapid method for the estimation of the glutamic-aspartic transaminase in tissues and its application to radiation sickness. *Arch Biochem* 1950; 28:36-42.
2. Reitman S and Frankel S. A colorimetric method for the determination of serum glutamic oxalacetic and glutamic pyruvic transaminases. *Am J Clin Pathol*, 1957; 28: 56-63.
3. Murray RL. Alanine aminotransferase. *In: Clinical Chemistry: Theory, Analysis, and Correlation*, 2<sup>nd</sup> ed. Kaplan LA, Pesce AJ, Eds. St. Louis: The C.V. Mosby Company, 1989: 895-898.
4. Wróblewski F, LaDue JS. Serum glutamic-pyruvic transaminase in cardiac and hepatic disease. *Proc Soc Exp Biol Med*, 1956; 91: 569-571.
5. Bergmeyer HU, Horder M. IFCC methods for the measurement of catalytic concentration of enzymes. Part 3. IFCC method for alanine aminotransferase. *J Clin Chem Clin Biochem*, 1980; 18: 521-534.
6. Howe PE. The use of sodium sulfate as the globulin precipitant in the determination of proteins in blood. *J Biol Chem*, 1921; 49: 93-107.
7. Howe PE. The determination of proteins in blood — a micro method. *J Biol Chem*, 1921; 49: 109-113.
8. Wolfson WQ, et al. A rapid procedure for the estimation of total protein, true albumin, total globulin, alpha globulin, beta globulin and gamma globulin in 10 ml of serum. *Am J Clin Pathol*, 1948; 18: 723-730.
9. Saifer A, Gerstenfeld S, Vacsler F. Photometric microdetermination of total serum globulins by means of a tryptophan reaction. *Clin Chem*, 1961; 7: 626-636.
10. Saifer A, Marven T. The photometric microdetermination of serum total globulins with a tryptophan reaction: a modified procedure. *Clin Chem*, 1966; 12: 414-417.
11. Gendler SM. Albumin. *In: Clinical chemistry: Theory, Analysis, and Correlation*. LA Kaplan and Pesce AJ, Eds. St. Louis: The C.V. Mosby Company, 1989: 1029-1033.
12. Webster D, Bignell AHC, Attwood EC. An assessment of the suitability of bromocresol green for the determination of serum albumin. *Clin Chim Acta* 53, 1974; 101-108.
13. Louderback A, Mealey EH, Taylor NA. A new dye-binding technique using bromocresol purple for determination of albumin in serum. *Clin Chem* 14, 1968; 793-794. (Abstract)
14. Pinnell AE, Northam BE. New automated dye-binding method for serum albumin determination with bromocresol purple. *Clin Chem*, 1978; 24: 80-86.
15. King EJ, Armstrong AR. A convenient method for determining serum and bile phosphatase activity. *Can Med Assoc J*, 1931; 31: 376-381.
16. Kind PRN, King EJ. Estimation of plasma phosphatase by determination of hydrolysed phenol with amino-antipyrine. *J Clin Pathol*, 1954; 7: 322-326.
17. Ohmori Y. Uber die Phosphomonoesterase. *Enzymologia*, 1937; 4: 217-231.
18. Fujita H. Umber die Mikrobestimmung der Blutphosphatase. *J Biochem, Japan*, 1939; 30: 69-87.

### 13. Bibliography (continued)

19. Petittlerc C et al. Mechanism of action of  $Mg^{2+}$  and  $Zn^{2+}$  on rat placental alkaline phosphatase. I. Studies on the soluble  $Zn^{2+}$  and  $Mg^{2+}$  alkaline phosphatase. *Can J Biochem*, 1975; 53: 1089-1100.
20. Tietz NW, et al. A reference method for measurement of alkaline phosphatase activity in human serum. *Clin Chem*, 1983; 29: 751-761.
21. Bowers GN, et al. IFCC methods for the measurement of catalytic concentration of enzymes. Part I. General considerations concerning the determination of the catalytic concentration of an enzyme in the blood serum or plasma of man. *Clin Chem Acta*, 1979; 98: 163F-174F.
22. Karmen A. A note on the spectrophotometric assay of glutamic-oxalacetic transaminase in human blood serum. *J Clin Invest*, 1955; 34: 131-133.
23. Bergmeyer et al. Provisional recommendations on IFCC methods for the measurement of catalytic concentrations of enzymes. Part 2. IFCC method for aspartate aminotransferase. *Clin Chem*, 1977; 23: 887-899.
24. Bergmeyer et al. Provisional recommendations on IFCC methods for the measurement of catalytic concentrations of enzymes. Part 2. Revised IFCC method for aspartate aminotransferase. *Clin Chem*, 1978; 24: 720-721.
25. Van Den Bergh AAH, Müller P. Über eine Direkte und eine indirekte Diazoreaktion auf Bilirubin. *Biochem Z*, 1916; 77: 90-103.
26. Malloy HT, Evelyn KA. The determination of bilirubin with the photoelectric colorimeter. *J Biol Chem*, 1937; 119: 481-490
27. Dumas BT, Wu T-W. The measurement of bilirubin fractions in serum. *Crit Rev Clin Lab Sci*, 1991; 28: 415-445.
28. Murao S and Tanaka N. A new enzyme "bilirubin oxidase" produced by *Myrothecium verrucaria* MT-1. *Agric Biol Chem*, 1981; 45: 2383-2384.
29. Osaki S, Anderson S. Enzymatic determination of bilirubin. *Clin Chem*, 1984; 30: 971. (Abstract)
30. Perry B, et al. Measurement of total bilirubin by use of bilirubin oxidase. *Clin Chem*, 1986; 32: 329-332.
31. Dumas BT, et al. Measurement of total bilirubin by use of bilirubin oxidase. *Clin Chem*, 1986; 32: 1349-1353.
32. Otsuji S et al. A new enzymatic approach for estimating total and direct bilirubin. *Clin Biochem*, 1988; 21: 81-110.
33. Meites S. Bilirubin, direct reacting and total, modified Malloy-Evelyn method. *In: Faulkner WR and Meites S, eds. Selected methods of clinical chemistry.* Washington, DC: American Association for Clinical Chemistry, 1982; 9: 119-124.
34. Koller A and Kaplan LA. Total serum protein. *In: Clinical chemistry: Theory, Analysis, and Correlation*, 2<sup>nd</sup> ed. Kaplan LA, Pesce AJ, Eds. St. Louis: The C.V. Mosby Company, 1989: 1057-1060.
35. Reigler E. Eine kolorimetrische Bestimmungsmethode des Eiweisses. *Z Anal Chem*, 1914; 53: 242-245.
36. Weichselbaum TE. An accurate and rapid method for the determination of proteins in small amounts of blood serum and plasma. *Am J Clin Path*, 1946; 16: 40-49.
37. Dumas BT et al. A candidate reference method for determination of total protein in serum. I. Development and validation. *Clin Chem*, 1981; 27: 1642-1650.
38. National Committee for Clinical Laboratory Standards. Physician's office laboratory guidelines, tentative guideline, 2<sup>nd</sup> ed. NCCLS document POL1-T2. Wayne, PA: NCCLS, 1992
39. National Committee for Clinical Laboratory Standards. Procedures for the handling and processing of blood specimens; approved guideline, 2<sup>nd</sup> ed. NCCLS document H18-T. Wayne, PA: NCCLS, 1988.
40. Rehak NN, Chiang BT. Storage of whole blood: effect of temperature on the measured concentration of analytes in serum. *Clin Chem*, 1988; 34: 2111-2114.
41. Sherwin JE, Obernolte R. Bilirubin. *In: Clinical chemistry: theory, analysis, and correlation*, 2<sup>nd</sup> ed. Kaplan LA and Pesce AJ, Eds. St. Louis: The C.V. Mosby Company, 1989:1009-1015.
42. Henry RJ, Cannon DC, Winkelman JW. *Clinical chemistry: Principles and techniques*, 2<sup>nd</sup> ed. New York: Harper and Row; 1974; 417- 421, 1058 -1059.
43. Benet LZ, Williams RL. Design and optimization of dosage regimens: pharmacokinetic data. *In: Gilman AG, et al. Eds. Goodman and Gilman's The Pharmacological Basis of Therapeutics*, 8<sup>th</sup> ed. New York: McGraw-Hill, Inc., 1990; 1650-1735.
44. Moss DW, Henderson AR. Enzymes. *In: Tietz Textbook of Clinical Chemistry*, 2nd ed. Burtis CA and Ashwood ER, Eds. Philadelphia. W.B. Saunders Company, 1994; 735-896.
45. National Committee for Clinical Laboratory Standards. Interference testing in clinical chemistry; approved guideline. NCCLS document EP7-A. Wayne, PA: NCCLS, 2002.
46. Painter PC, Cope JY, Smith JL. Reference Information for the clinical laboratory. *In: Tietz Textbook of Clinical Chemistry*, 3<sup>rd</sup> ed. Burtis CA, Ashwood ER, Eds. Philadelphia. WB Saunders Company, 1999: 1788-1846.
47. Painter PC, Cope JY, Smith JL. Appendix. *In: Tietz Textbook of Clinical Chemistry*, 2<sup>nd</sup> ed. Burtis CA, Ashwood ER, Eds. Philadelphia. W.B. Saunders Company, 1994: 2161-2217.
48. Young DS. Effects of drugs on clinical laboratory tests, 3<sup>rd</sup> ed. Washington, DC: AACC Press, 1990.
49. National Committee for Clinical Laboratory Standards. Evaluation of precision performance of clinical chemistry devices, 2nd ed.; approved guideline. NCCLS document EP5-A. Wayne, PA: NCCLS, 1999.