



piccolo^{xpress}[®] chemistry analyzer

FAST, EASY, ACCURATE ON-SITE CHEMISTRY TESTING

REAL-TIME CLIA-WAIVED CHEMISTRY

EXCEPTIONALLY VERSATILE

Wide range of **CLIA waived** panels to increase your on-site diagnostic capabilities

- **Comprehensive Metabolic Panel**
ALB, ALP, ALT, AST, BUN, Ca, Cl⁻, CRE, GLU, K⁺, Na⁺, TBIL, tCO₂, TP
- **Basic Metabolic Panel**
BUN, Ca, Cl⁻, CRE, GLU, K⁺, Na⁺, tCO₂
- **Lipid Panel**
CHOL, CHOL/HDL*, nHDLc*, HDL, LDL*, TRIG, VLDL*
- **Lipid Panel Plus**
ALT, AST, CHOL, CHOL/HDL*, nHDLc*, GLU, HDL, LDL*, TRIG, VLDL*
- **Liver Panel Plus**
ALB, ALP, ALT, AMY, AST, GGT, TBIL, TP
- **General Chemistry 6**
ALT, AST, BUN, CRE, GGT, GLU
- **General Chemistry 13**
ALB, ALP, ALT, AMY, AST, BUN, Ca, CRE, GGT, GLU, TBIL, TP, UA
- **Electrolyte Panel** - Cl⁻, K⁺, Na⁺, tCO₂
- **Kidney Check** - BUN, CRE
- **Renal Function Panel**
ALB, BUN, Ca, Cl⁻, CRE, GLU, K⁺, Na⁺, PHOS, tCO₂
- **MetLyte 8 Panel**
BUN, CK, Cl⁻, CRE, GLU, K⁺, Na⁺, tCO₂

- **Hepatic Function Panel**
ALB, ALP, ALT, AST, DBIL, TBIL, TP
- **Basic Metabolic Panel Plus**
BUN, Ca, Cl⁻, CRE, GLU, K⁺, Na⁺, tCO₂, Mg, Lactate Dehydrogenase
- **MetLyte Plus CRP**
BUN, CK, Cl⁻, CRE, GLU, K⁺, Na⁺, tCO₂, CRP
- **BioChemistry Panel Plus**
ALB, ALP, ALT, AMY, AST, BUN, Ca, CRE, CRP, GGT, GLU, TP, UA
- **MetLac 12 Panel**
ALB, BUN, Ca, Cl⁻, CRE, GLU, LAC, K⁺, Mg, Na⁺, Phos, tCO₂

EXCEPTIONALLY FAST

From sample to complete panel results in 3 simple steps and approximately 12 minutes



1 / Add sample



2 / Insert disc



3 / Read results

EXCEPTIONALLY EASY

Fully automated system requires no special operating skills

- 0.1cc whole blood, serum or plasma
- Self calibrates with every test run - no manual data entry
- Full-color touch screen with intuitive user interface
- Results are printed on a self-adhesive hard copy or via electronic transfer, allowing immediate access to critical information at the point-of-care
- Large memory capacity stores up to 5000 patient records
- Externally accessible bi-directional USB data transfer technology facilitates EMR connectivity

EXCEPTIONALLY ACCURATE

Advanced technology ensures sound results

- Centrifugal action mixes the reagents and sample in a disposable self-contained disc – reducing procedural errors
- Accuracy and precision comparable to larger, more expensive laboratory analyzers
- iQC[®] system performs continuous quality control checks on every run, eliminating daily control testing requirements
- Every self-contained reagent disc is barcoded for automatic self-calibration


CLIA WAIVED

*Calculated



3240 Whipple Rd., Union City, CA 94587 | 800.822.2947 | www.piccoloxpress.com

©2015 Abaxis, Inc., Piccolo Xpress[®] and iQC[®] are registered trademarks of Abaxis, Inc. 888-3010 Rev.K

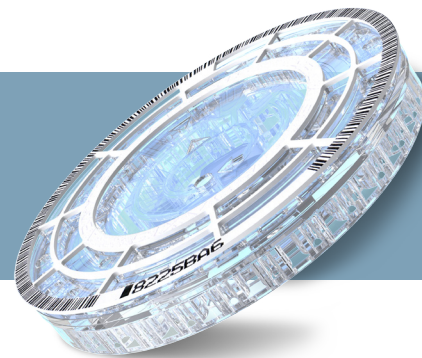

PICCOLO XPRESS PANEL DETAIL

Analyte	CLIA WAIVED												MODERATELY COMPLEX			
	Comprehensive Metabolic Panel	Basic Metabolic Panel	Liver Panel Plus	General Chemistry 6	General Chemistry 13	Kidney Check	Renal Function Panel	Lipid Panel	Lipid Panel Plus	Electrolyte Panel	MetLyte 8 Panel	Hepatic Function Panel	Basic Metabolic Panel Plus	MetLyte Plus CRP	Biochemistry Panel Plus	MetLac 12 Panel
ALB	●		●		●		●					●			●	●
ALP	●		●		●							●			●	
ALT	●		●	●	●			●							●	
AMY			●		●										●	
AST	●		●	●	●			●				●			●	
BUN	●	●		●	●	●	●			●		●	●	●	●	●
Ca	●	●			●		●					●			●	●
CK										●			●			
Cl ⁻	●	●					●		●	●		●	●			●
CRE	●	●		●	●	●	●			●		●	●	●	●	●
GGT			●	●	●										●	
GLU	●	●		●			●		●		●	●	●	●	●	●
K ⁺	●	●					●		●	●		●	●			●
Na ⁺	●	●					●		●	●		●	●			●
PHOS							●									●
TBIL	●		●		●							●				
tCO ₂	●	●					●		●	●		●	●			●
TP	●		●		●							●			●	
UA					●										●	
CHOL								●	●							
CHOL/HDL*								●	●							
HDL								●	●							
LDL*								●	●							
nHDLc*								●	●							
TRIG								●	●							
VLDL*								●	●							
CRP													●	●		
DBIL												●				
LAC																●
LDH												●				
Mg												●				●

* Calculated